


Realize Language: What do the codes mean in the Log page Log column?

When someone uses an AAC device or software, it's possible to log what happens. This might include pressing a button to speak a word, clearing the display, switching to a new page, turning down the volume, and any other action. We call all these *events* and the basic idea behind data logging as a feature is that it is a way to track and record events.

And usually when an event takes place the AAC software creates a special code to mark what happened. These codes may vary slightly from device to device but the Realize Language software can recognize these, which is why it works with a range of software and devices.

Here are list of codes for the specific devices and software. In order to see these you need to be on the **Log** page and set for the *Full Layout* with Details showing:

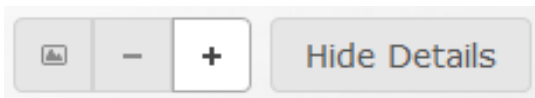

Note that if you hit the "Hide Details" button you will hide some of the codes, so if you want to explore them all, you need to make sure you are seeing "Hide Details" just as in the image above.

Prentke Romich devices

CODE	MEANING	EXAMPLE		
ACT	A selection was made from an ACTivity Row.	ACT	17:47:2	<i>circle</i>
CTL	A non-language key activation.	CTL	14:56:25	Page transition from QWERTY to MESSAGES.
DWP	Derived Word Prediction; a word generated by being predicted after a regular WPR word prediction.	SPE WPR DPR	11:31:33 11:31:35 11:31:37	<i>s</i> <i>speak</i> <i>speakers</i>
LOC	LOCation of a key by Row + Column and [Name of key].	LOC LOC ACT	15:59:46 15:59:48 15:59:50	I5 [MORNING] E1 [WATER BODIES] B1 <i>beach</i>
PAG	A selection was made from a PAGE key.	CTL PAG PAG	15:09:46 15:09:51 15:09:56	Page transition from SPELL to NUMBERS C2 <i>one</i> E5 <i>forty</i>
SEM	A key was used in a SEMantic Compaction <i>sequence</i> .	LOC SEM LOC SEM LOC LOC SEM	15:58:46 15:58:47 15:58:49 15:58:50 15:58:52 15:58:54 15:58:56	E5 [DOG] C3 [VERB] <i>come</i> G2 [CONJ] G2 [CONJ] <i>and</i> K6 [TRUCK] E7 [FLY] C3 [VERB] <i>fly</i>
SPE	Character is generated by SPELLing a single letter/number.	SPE SPE SPE SPE SPE	09:39:56 09:39:57 09:39:58 09:39:59 09:40:04	D3 [] <i>s</i> I2 [] <i>u</i> F2 [] <i>r</i> E2 [] <i>e</i> A2 [SPACE]
WPR	Word Prediction:	SPE SPE WPR	10:26:35 10:26:36 10:26:38	C5 <i>s</i> H4 <i>u</i> F3 <i>supposed</i>

TouchChat/LAMP Words for Life Software

CODE	MEANING	EXAMPLE		
CTL	A non-language key activation. Not visible when "Hide details" is selected on the Realize Log page.	CTL	14:55:45	EnableLAM: ON
		CTL	14:56:25	Vocab: Copy LAMP Full
		CTL	14:58:30	PGG (Nav Home): HOME
DEL	Indicates something has been DELETED, such as a character, word, or contents of the display.	LOC	17:01:58	0,2,1
		DEL	17:01:58	DISPLAY
LOC	LOCation of a key by Type (typically a 0), Row, and Column. Not visible when "Hide details" is selected on the Realize Log page.	LOC	15:59:46	0,2,1
		LOC	15:59:48	0,4,1
PAG	A selection was made from a PAGE key.	PAG	15:09:51	0,2,2 <i>help</i>
PGG	Follows a CTL code and specifies the type of Navigation that is happening e.g. to Home, Page, or temporary Visit.	CTL	19:25:54	PGG (Nav Home): HOME
		CTL	19:25:58	PGG (Visit): Groups
		CTL	19:26:03	PGG (Nav Back):
SPE	Character is generated by SPELLing a single letter/number.	SPE	11:31:33	0,7,7 <i>c</i>
		SPE	11:31:35	0,5,1 <i>a</i>
		SPE	11:31:37	0,3,7 <i>t</i>
SPK	Indicates what is actually spoken by the device, which may be different from the text on a button.	PAG	17:00:53	0,5,4 <i>it</i>
		SPK	17:00:53	<i>it</i>
		CTL	17:00:54	PGG (Visit): IT WORDS
		PAG	17:00:54	0,2,3 <i>will</i>
		SPK	17:00:54	<i>will</i>

Salttillo NovaChat Devices

CODE	MEANING	EXAMPLE
ACTION	Specifies a non-speech ACTION that happens when a key is selected. Not visible when "Hide details" is selected on the Realize Log page.	PAG 15:09:51 / CTL 15:09:56 ACTION (Speech Message) CTL 15:09:51 ACTION <Visit 'Emotion Page' PAG 15:09:56 <i>feel</i> CTL 15:09:56 ACTION (Speech Message) PAG 15:10:02 <i>happy</i>
CTL	A non-language key activation. Not visible when "Hide details" is selected on the Realize Log page.	CTL 14:56:25 ACTION (Static Word Prediction,2)
LOC	LOCation of a key by Row + Column. Not visible when "Hide details" is selected on the Realize Log page.	LOC 15:59:46 2,1 LOC 15:59:48 4,1
PAG	A selection was made from a PAGE key.	PAG 15:09:51 2,2 <i>help</i>
SPE	Character is generated by SPELLing a single letter/number.	SPE 14:16:05 5,6 <i>h</i> CTL 14:16:05 ACTION(Speech Message) SPE 14:16:07 3,9 <i>o</i> CTL 14:16:07 ACTION(Speech Message) SPE 14:16:10 3,2 <i>w</i> CTL 14:16:10 ACTION(Speech Message) SPE 14:16:12 9,1 CTL 14:16:12 ACTION(Speech Message)